

Miles-Made

Pentex® Bovine Albumins...Exactly What You Need.

Don't settle for something "close"—Miles probably has exactly what you need. Choose from the most comprehensive selection of forms, sizes, and derivatives of bovine albumin available. We offer 18 different catalog varieties, both powder and solution, in convenient packages from grams to kilograms. Look at this list...it's exactly what you've been looking for.

New! Sulfhydryl Modified—chemically modified to reduce free sulfhydryl groups to <1%. (Code No. 81-018)

Monomer (powder)—>99% monomer when manufactured. (Code No. 81-028)

Crystallized (powder)—the high-purity albumin for stringent research systems. (Code No. 81-001)

Fraction V, Fatty Acid Free (powder)—has the largest number of available binding sites. (Code No. 82-002)

Fraction V (powder)—standard general-purpose albumin. (Code No. 81-003)

Polymer Enhanced, 30% Solution—naturally polymerized, high avidity for serology. (Code No. 81-078)

Albumin: An Overview and Bibliography—an up-to-date monograph, featuring the most comprehensive bibliography (over 800 references) on the subject. For your copy and to order Pentex Bovine Albumins, call or write:

Miles Scientific
Div. Miles Laboratories, Inc.
30W475 North Aurora Road
Naperville, IL 60566
USA
Tel. Toll-free 800-348-7465
(in Illinois, 312-357-3720,
in Canada, 416-248-0771)

Miles Scientific
Div. Miles Laboratories, Ltd.
Stoke Court, Stoke Poges
Slough SL2 4LY
United Kingdom
Tel. Farnham Common 5151



1981 MILES LABORATORIES, INC.

MRP 5146812

State-of-the-art technology in clinical microbiology

RAPID METHODS AND AUTOMATION IN MICROBIOLOGY

Editor: Richard C. Tilton

This exciting new book, based on papers presented at the Third International Symposium on Rapid Methods and Automation in Microbiology, May 1981, deals with the following areas:

- Rapid Methods and Automation in Microbiology—1981
- Microbiological Approaches to the Rapid Detection of Microorganisms in Body Fluids
- Hydrophobic Grid-Membrane Microbiology—Principles and Practice
- A Decade of Automation
- Concept of Rapid Reporting in Microbiology
- Immunological Diagnosis of Infectious Diseases
- Chromatographic Analysis: a New Future for Clinical Microbiology
- Microbiological Approaches to Rapid Detection of Microorganisms in Urine
- Rapid Methods in Food Microbiology
- Systems Approaches to Diagnostic Microbiology
- Rapid Identification of Viral, Chlamydial, and Mycoplasmic Agents
- Current Methods for Rapid Identification of Clinically Important Anaerobes
- Applications of the Autobac System
- Clinical and Research Applications of the MS-2 System
- Nonautomated Rapid Methods in Microbiology

ASM

**American Society
for Microbiology**

ORDERING INFORMATION:

Publication Date: February 1982. 361 pages. Clothbound. Author and subject indexes. ISBN: 0-914826-38-7. \$33.00. Price subject to change without notice.

Order from the publisher:

Publication Sales, 1913 I Street, NW, Washington, DC 20006

Meeting the bacteriologist's need for a compact handbook of laboratory methods—

Manual of Methods for General Bacteriology


EDITOR-IN-CHIEF:
Philipp Gerhardt

SECTION EDITORS:
Ralph N. Costilow
Noel R. Krieg
R. G. E. Murray
Eugene W. Nester
G. Briggs Phillips
Willis A. Wood

Published and distributed by:



American Society for Microbiology
1913 I Street, NW
Washington, DC 20006

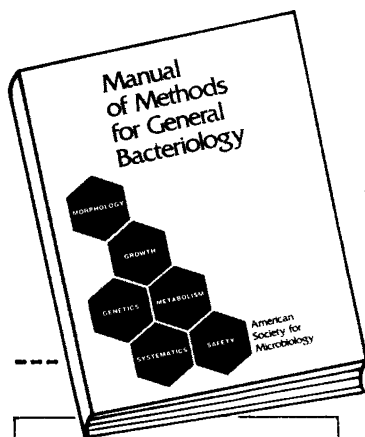

This *Manual* covers all the taxa of bacteria, and its approach is generalized for both pure and applied bacteriology. Aimed at senior undergraduate and beginning graduate science majors, at teachers and researchers, the *Manual* provides practical descriptions of laboratory procedures that are known to be reliable and that are in frequent use. It is an invaluable tool for all involved in bacteriology in the laboratory setting.

The book comprises six sections, each of which is self-contained, with its own table of contents to help the reader find specific methods as easily as possible. The contents are listed below.

METHODOLOGY FOR GENERAL BACTERIOLOGY

MORPHOLOGY—Light Microscopy—Specimen Preparation for Light Microscopy—Determinative Methods of Light Microscopy—Electron Microscopy—Cell Fractionation

GROWTH—Biophysical Factors in Growth—Biochemical Factors in Growth—Enrichment and Isolation—Solid Culture—Liquid Culture—Growth Measurement—Preservation

GENETICS—Gene Mutation—Gene Transfer—Plasmids

METABOLISM—Physical Methods—Chemical Composition—Enzymatic Activity—Permeability and Transport

SYSTEMATICS—General Characterization—Nomenclature—Taxonomy—Genetic Characterization

LABORATORY SAFETY—Sterilization—Contamination and Disinfection

ADDENDUM—Diluents and Biomass Measurement

ORDERING INFORMATION

Publication date: March 1981 525 pages, fully illustrated

Clothbound

Flexible binding

ISBN: 0-914826-29-8

ISBN: 0-914826-30-7

\$30.00 (\$22.00 to members of ASM) \$26.00 (\$19.00 to members of ASM)

Prices subject to change without notice

Please send me _____ copy(ies) of
MANUAL OF METHODS FOR GENERAL BACTERIOLOGY
for \$30.00 (Cloth); _____ \$26.00 (Flexible)

Send book(s) to: _____ I have enclosed \$ _____

Name _____

Address _____

City/State/Zip _____

I am a member of ASM.
Send me _____ copy(ies):

_____ (limit of 3 at member prices)

\$22.00 (Cloth)

\$19.00 (Flexible)

I have enclosed \$ _____

The newest addition to the annual
MICROBIOLOGY series,

MICROBIOLOGY-

1982 puts into focus
exciting new work in
the following major areas of microbiological research:

CELL-CELL INTERACTIONS IN PATHOGENESIS

Bacterial Adhesion in Pathogenesis

- Introduction
- Autochthonous Populations
- Role of Surface Proteins Involved in Surface Adhesion
- Role of Surface Polysaccharides in Bacterial Adhesion
- Variable Surface Factors in Bacterial Adhesion

The Macrophage in Host Defense

Cell Surface Receptors Involved in the Immune System

COMPARATIVE MOLECULAR BIOLOGY

Streptococcal Genetics

- Conjugation
- Transformation, Transduction, and Bacteriophage
- Drug Resistance
- Genetics and Physiology
- Metabolic and Bacteriocin-Related Plasmids
- Development and Use of Recombinant DNA Technology

Bacillus subtilis Chromosome

- Structure and Modification
- Molecular Cloning

As usual in this well-established series, these extended abstracts are carefully selected for their importance to the development of the knowledge base in each of these areas. The book will be particularly useful in research laboratories and institutional libraries.

Bibliographic and Ordering Information

Editor: David Schlessinger.

July 1982, 420 pages, clothbound. \$15.00 to ASM members; \$25.00 to others. Please note that ASM members must remit payment with order. Orders from outside the United States and Canada must be prepaid before shipment. Prices are subject to change without notice. ISBN: 0-914826-42-5. Standing orders to the series are welcome.

Previous volumes in the series, which started publication in 1975 with the volume entitled *Microbiology — 1974*, are still available in limited quantities. Prices for each preceding volume are \$15.00 for ASM members and \$25.00 for others.

ASM
American Society
for Microbiology

Published and distributed by:
American Society for Microbiology
Publications Office
1913 I Street, NW
Washington, DC 20006

An enlightening discussion of

HOST-PARASITE INTERACTIONS IN PERIODONTAL DISEASES

EDITORS: Robert J. Genco and Stephan E. Mergenhagen

The papers in this book were first presented at a symposium held in Buffalo, New York, in May 1981. They are divided into six sections:

- I. Etiology and Microbial Specificity
- II. Colonization by Periodontopathic Organisms
- III. Pathogenic and Virulence Factors in Periodontal Diseases
- IV. Cellular Aspects of Host Responses
- V. Effector Systems—Specific and Nonspecific
- VI. Clinical Implications of Host-Parasite Interactions

ORDERING INFORMATION:

Publication date: January
1982. 414 pages. Clothbound.
Author and subject indexes.
ISBN: 0-914826-37-9. \$27.00
Price subject to change
without notice.

Order from the publisher:



**AMERICAN SOCIETY
FOR MICROBIOLOGY**

Publication Sales
1913 I Street, NW
Washington, DC 20006

The world's two most widely used
clinical laboratory manuals are now
available in completely revised
and updated new editions:

Manual of Clinical Microbiology,

3RD EDITION
JULY 1980

Editors: Edwin H. Lennette, Albert Balows, William J. Hausler, Jr., Joseph P. Truant
1,047 pages, fully illustrated.
Clothbound, \$30.00 (\$22.00 to ASM Members)
Flexible binding, \$26.00 (\$19.00 to ASM Members)

Manual of Clinical Immunology,

2ND EDITION
JULY 1980

Editors: Noel R. Rose and Herman Friedman
1,105 pages, fully illustrated.
Clothbound, \$30.00 (\$22.00 to ASM Members)
Flexible binding, \$26.00 (\$19.00 to ASM Members)

Order your copies today from the publisher:
American Society for Microbiology
1913 I Street, NW • Washington, DC 20006



ASM
American Society
for Microbiology



RIBI
IMMUNOCHEM
RESEARCH, INC.

REAGENTS
FOR BASIC
IMMUNOLOGY
RESEARCH

STM *B-lymphocyte* MITOGEN

S. typhimurium Mitogen (STM) is a proteinacious B-cell mitogen prepared from purified cell walls of the heptosless Re-mutant of *S. typhimurium*. STM is non-toxic as measured by chick embryo lethality (CELD₅₀ > 10 μ g) and lymphocytes *in vitro* (dosages of 1 mg/ml). STM has extremely high mitogenic activity for B-lymphocytes from both endotoxin responder and non-responder mice, nude mice and rats.

Mitogenic Activity of STM, PHA and *E. coli* LPS on spleen cells.*

Cultured cells Tested	³ H-thymidine incorporation (CPM)			
	STM 50 μ g/ml	PHA 1 μ g/ml	LPS 50 μ g/ml	None
C57BL/10 ScN Endotoxin non-responder	18,588	20,288	1,109	695
(C57BL/6 \times DBA/2)F ₁ Endotoxin responder	56,998	43,667	11,766	769
BALB/c nu/+	70,544	53,940	11,440	1,030
BALB/c nu/nu	36,283	902	11,719	751
Rat B-cells	57,600 (10 μ g/ml)	—	25,700 (10 μ g/ml)	1,700
Rat T-cells	2,107 (10 μ g/ml)	—	1,179 (10 μ g/ml)	696

* 1 \times 10⁶ cells/ml of each strain were cultured with or without mitogens in 0.2ml of RPMI-1640 (5% FCS) in 96 well microtiter plates. One μ ci ³H-thymidine was added to each well the last 18 hours of a 48 hour incubation and ³H-incorporation was measured by standard scintillation techniques. Results of triplicate samples from each of three mice per group are expressed as the mean CPM.

STM is your choice for studying the mitogenic activity and polyclonal activation of B-lymphocytes, and B-cell suppression *in vitro*.

FOR FURTHER INFORMATION,
RIBI IMMUNOCHEM RESEARCH, INC.

P.O. BOX 1409, HAMILTON, MICHIGAN 49840 USA
TELEPHONE (406) 363-6214